

HERBST


Bite jumping hinges.
Safe. Proven. Effective.

An overview of our Herbst bite jumping hinges.




The classic correction of distal occlusions.

Dentaurum's program consists of 4 versions of the proven Herbst bite jumping hinge, whereby each version has an excellent finish. Even at high stress levels, Dentaurum guarantees that there are sufficient safety reserves to ensure therapy remains effective. Since the appliance is fixed directly on the teeth, a force is constantly exerted without having to involve the patient. Results can be seen after only a short time.

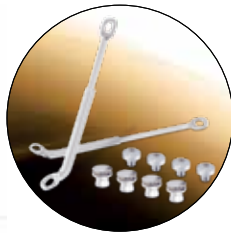


Herbst I - the classic version.

The Herbst I bite jumping hinge is the classic version that is welded onto bands or cast splints, depending on the aim of the treatment and the experience of the orthodontist.

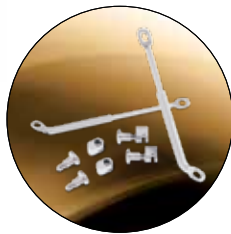

Herbst II - the modified version.

Here, the classic bite jumping hinge has been modified so that it can be quickly and easily screwed onto the arches of multi-bracket appliances. This way, it can be coupled with the appliance already in use.

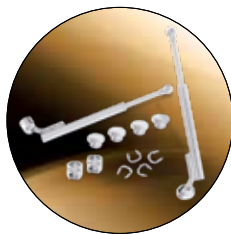

Herbst IV - the comfortable version.

The ball and socket joint offers the patient greater lateral movement, making the appliance noticeably more comfortable. The bite jumping hinge is fixed to the base with c-clips, thereby avoiding the use of screws.



Herbst TS - the compact version.

This model is fitted with an internal telescope, preventing the hinge from sliding apart (e.g. if patient opens his/her mouth wide). Its compact design makes Herbst TS the best solution for small jaws. The new laser markings on the telescopic element indicating size and side guarantee clear allocation to the jaw.

The advantages at a glance.



Herbst I set (REF 607-120-00)
on bands



F-bases (concave shape) for
Herbst I / Herbst TS
(REF 607-127-30)

Hexagon socket screws for
Herbst I / Herbst TS
(REF 607-122-30)



Hex key for Herbst I / II / TS for easier chairside
treatment (REF 607-144-00)

Fast therapy.

As a rule, it takes 6 to 8 months to achieve a neutral relation. The best time to start treatment is during the growth spurt in adolescence. The subsequent retention period with a removable appliance will then also be short.

Fast. Reliable. Time-saving.

HERBST II

Spacer rings for Herbst I / II / IV
(lengths 1 – 5 mm)
(REF 607-103-00, 1 mm)



Hexagon socket screws with long shank (REF 607-108-20) and nut (REF 607-111-00) for Herbst II



Herbst II set (REF 607-108-00)
on remanium® arch, rectangular

Reliable therapy.

For patients with poor compliance, this fixed appliance offers the advantage that the therapy is effective around the clock even without patient cooperation.

Time-saving combined therapy.

Combining Herbst bite jumping hinges with multi-bracket appliances allows different treatments to be carried out parallel to one another, saving time and also money.

A program for the highest demands.

HERBST IV



Herbst IV set (REF 607-115-00)
on bands



Telescopic tube with ball and
socket joint (REF 607-125-00)

Bases for Herbst IV
(REF 607-118-00)



Lock washers (c-clips) for
Herbst IV (REF 607-116-00)

Exact therapy.

Since preparation is done on a model in the laboratory, it is possible to set the target occlusion optimally and with great precision. This is not only advantageous for the orthodontist, but also delivers results that are unrivaled in precision.

Effective therapy.

Herbst bite jumping hinges can also be used when the jaw position at the beginning of treatment is asymmetric. Treatment is therefore not only faster, but also more effective. This also increases the spectrum for which the Herbst appliances can be used.

Exact. Effective. Complete.

HERBST TS

Hexagon socket screws and bases for
Herbst I / Herbst TS
(REF 607-122-30 without base)



F-bases (concave shape) for
Herbst I / Herbst TS
(REF 607-127-30)



Herbst TS set, 22 mm
(REF 607-126-22) on bands



Spacer rings, crimpable
for Herbst TS (length 1 – 4 mm)
(REF 607-103-11, 1 mm)

Spacer rings (crimpable).

Spacer rings in different lengths are available for all Dentaaurum Herbst bite jumping hinges. They can therefore be put into use quickly, simply and safely when needed.

Complete program.

Accessories are available for each version of the Dentaaurum Herbst bite jumping hinges. Individual appliances can be prepared easily and in a timely manner using the accessories, and the right tool is always at hand.

Dentaurum

Germany | Benelux | España | France | Italia | Switzerland | Australia | Canada | USA
and in more than 130 countries worldwide.



DENTAURUM
QUALITY
WORLDWIDE
UNIQUE

	www.dentaurum.com		
			

Date of information: 2023-10
Subject to modifications



Turnstr. 31 | 75228 Ispringen | Germany | Phone +49 72 31/803-0 | Fax +49 72 31/803-295
www.dentaurum.com | info@dentaurum.com