

EN



OUR
3D RETAINER
IS A PERFECT FIT



D
DENTAURUM
1886



THE

3D-RETAINER

AS INDIVIDUAL
AS EVERY PATIENT



DEVELOPED BY
EXPERTS

The **prime4me[®] RETAIN3R** is
the tailored solution for your patients

Thanks to the innovative 3D design, the retainer can be adapted to the **individual tooth structure**, even where space is tight.

Made of high-quality **titanium grade 5**, our 3D retainer offers not only a high level of biocompatibility, but is also ideal for patients who prefer nickel-free materials.

Since **the fit is so precise**, the retainer is comfortable to wear and easy to clean – making for happy patients and optimum maintenance of orthodontic treatment results. A life long!



More information

THE


3RD DIMENSION

MAKES THE
DIFFERENCE

**Flat, half-round
cross-sectional profile**


Delicate 3D retainer,
robust and durable,
thickness 0.5 mm

PRECISE · EXACT FIT · INDIVIDUAL · STABLE · DURABLE

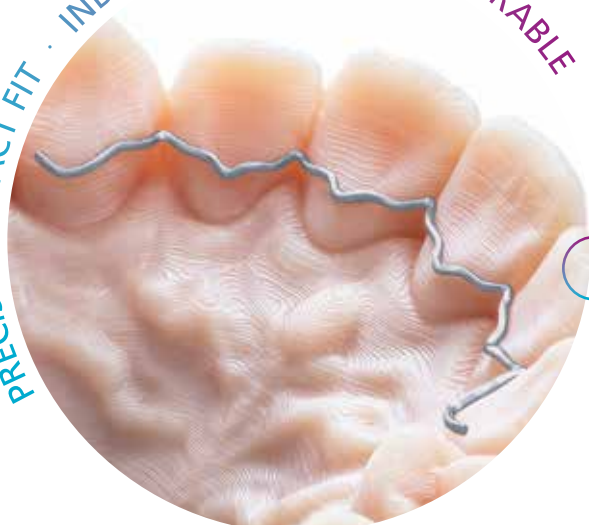

**Precision in
every dimension**

Sagittal, transverse and
vertical planes considered

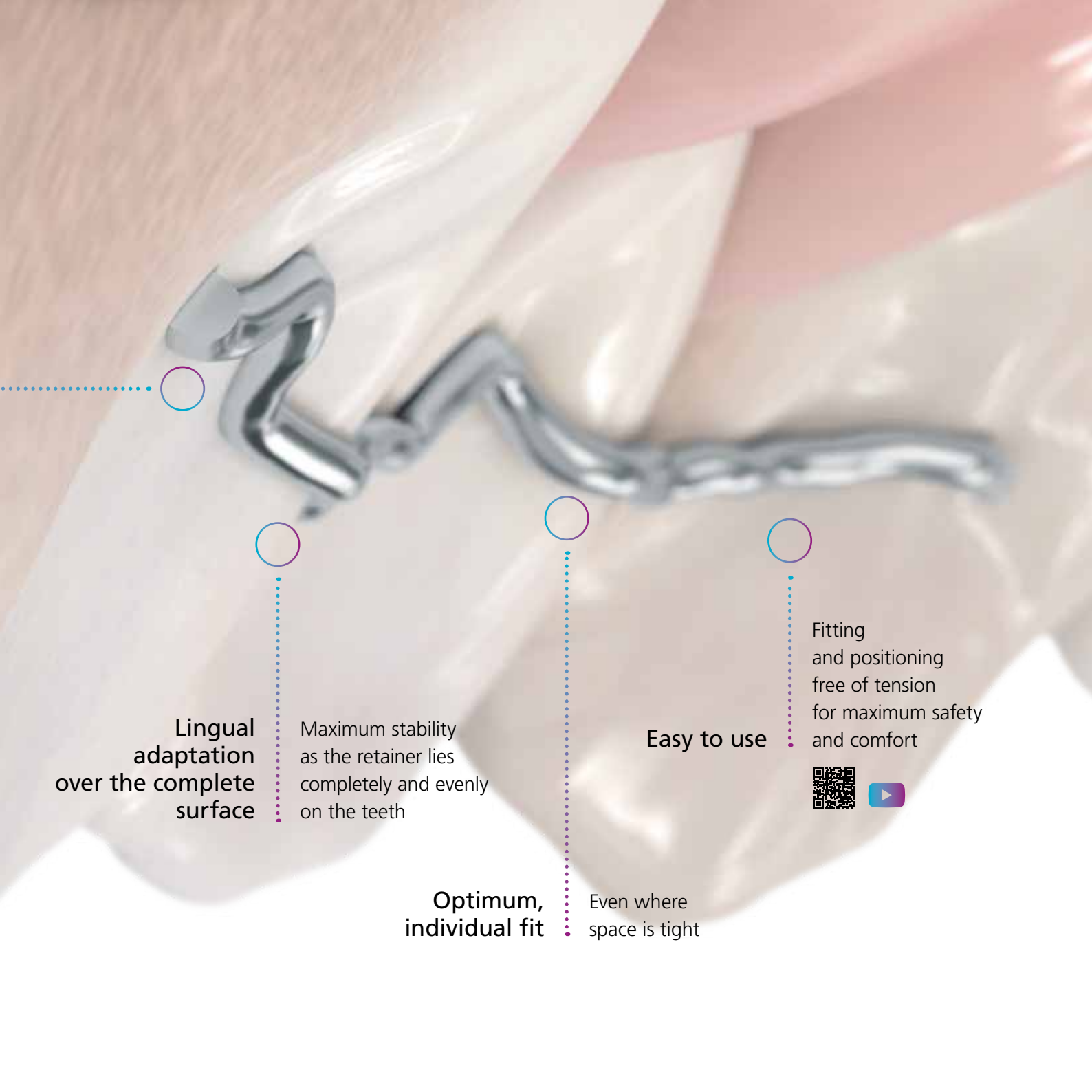

Lingual adaptation over the complete surface

Maximum stability as the retainer lies completely and evenly on the teeth

Optimum, individual fit

Even where space is tight

Easy to use

Fitting and positioning free of tension for maximum safety and comfort



Ti


Titanium

by DENTAURUM

High mechanical strength: ideal material for retainers.

High biocompatibility: optimum compatibility with the human body, also for patients who prefer nickel-free materials.

High corrosion resistance: improves product durability.

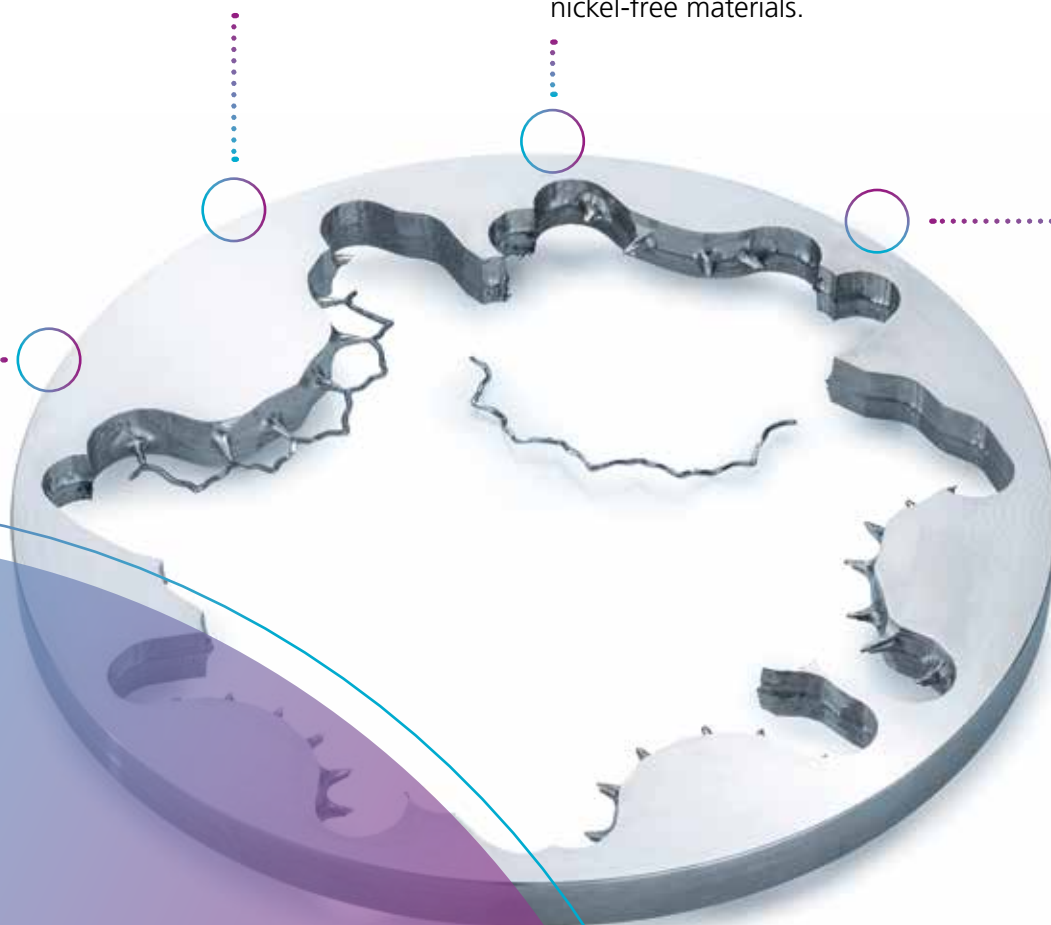

WE USE

TITANIUM GRADE 5

BENEFIT FROM OUR **KNOW-HOW**

From one piece:

no forming and
no bending.




Dentaurum high-tech

Our **prime4me® RETAIN3R** is passive thanks to a state-of-the-art milling technique: the material is machined and neither formed, nor bent. This production process is ideal for maintaining and stabilizing the results of your orthodontic treatment. The teeth still maintain their natural mobility. The bonding spots are sandblasted to improve the adhesive bond.

We offer 24 months guarantee.



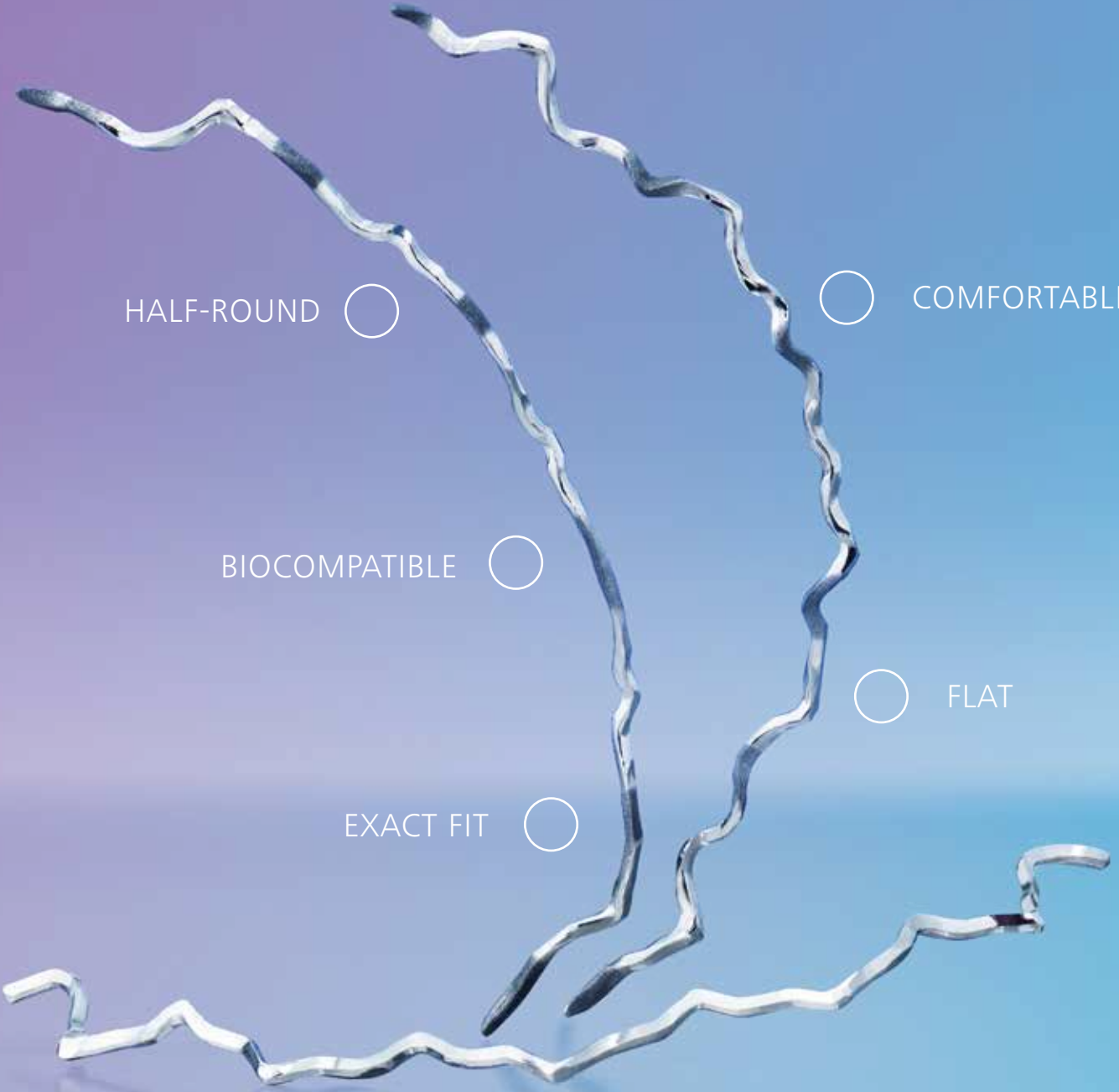

HALF-ROUND ○

COMFORTABLE ○

BIOCOMPATIBLE ○

FLAT ○

EXACT FIT ○

Comfortable to wear


The lingual retainer sits snugly, the bonding surfaces are therefore flatter, making it barely noticeable and comfortable to wear.

Improved oral hygiene

Oral hygiene is simple and effective thanks to the flat, half-round profile. This facilitates cleaning the teeth and the interdental spaces.

Highly compatible

Made of high-quality titanium grade 5, the retainer is also ideal for patients who prefer nickel-free materials.

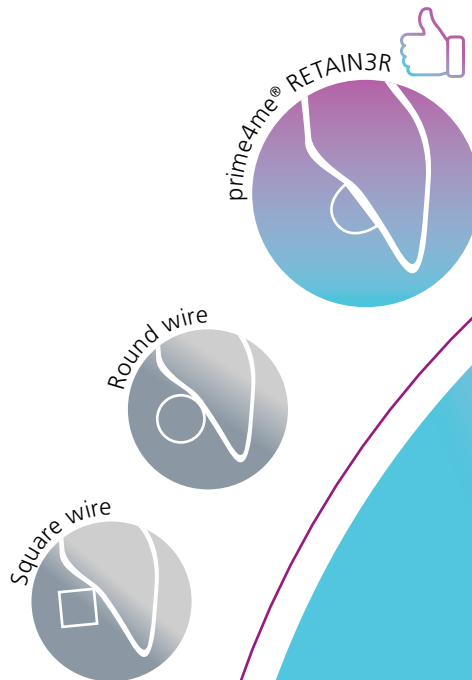

MANY
ADVANTAGES
FOR YOUR
PATIENTS

SIMPLE

TO ORDER

WITH A FEW
CLICKS



Here you can
order directly:
digital.dentaurum.de



1. Order with intra-oral
scan oder model scan



2. Approval of
optimum 3D design

Any questions?

Benefit from our CAD/CAM experts' competent advice.
We accompany you throughout the entire process.


Customer Support Digital:

☎ +49 72 31/803-280 🌐 digital@dentaurum.de

You can reach us (CET):

Monday thru Thursday: 8:00 a.m. to 5:00 p.m.

Friday: 8:00 a.m. to 3:00 p.m.

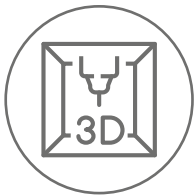

3. Production of
prime4me® RETAIN3R



4. Quality control
and cleaning



5. Handover to the
shipping department

.....
Delivery times depend on the country of
destination. You can find an exact overview here:



DURATION: **5 DAYS**



CONTACT
DENTAURUM

Dentaurum GmbH & Co. KG



Turnstr. 31
75228 Ispringen · Germany



info@dentaurum.com
www.dentaurum.com



+49 72 31/803-0



ONLINE SHOP
SHOP.DENTAURUM.COM